CA05FS

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© Sunnysoft s.r.o., distributor

DEVICE DESCRIPTION

The device is designed to detect the presence of alcohol in the human body. Before first use, carefully read this manual. Keep the manual for later reference.


DO NOT DRIVE UNDER THE INFLUENCE OF ALCOHOL!

 © Sunnysoft s.r.o., distributor

Display

The test result is shown in BrAC/per mille values corresponding to the 3rd decimal place in BAC.

The unit of measurement is set according to the given region.

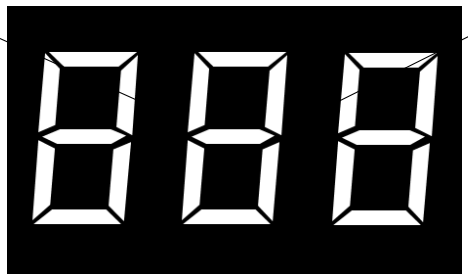The green LED indicates that the device is ready for use. The green LED also lights up when the measurement result is lower than the warning level.

READY

BATT

WARN

If the result is higher than 0.25 mg/l / 0.50 ‰ per mille, the red WARN LED lights up.

The yellow LED flashes when the battery is low. Replace the 2x 1.5 V alkaline batteries with new ones.

© Sunnysoft s.r.o., distributor

Operation

1. Insert the mouthpiece into the mouthpiece inlet.
2. Press the power on/off button to turn the device on.
 - The device first shows “Full display” and “on” on the display, while a welcome tone sounds. Then a countdown starts that lasts approximately 15 ∼ 35 s (warming up and preparing the device for use).
3. As soon as you hear a beep and “bLo” lights up on the display, start blowing into the mouthpiece until you hear a beep again. Then stop blowing into the device.
 - If you do not blow into the device within 20 s, the device turns off automatically. Turn the device on again.
4. For 3 ∼ 22 s a “bars” symbol will rotate, then the test result is shown for 20 s (mg/l BrAC or ‰ per mille) and an audible signal sounds.
 - If the test result is higher than 0.25 mg/l / 0.50 ‰ (default setting), the “WARN” signal (red LED) lights up and a warning “Alarm” sounds.
 - The measurement result is shown on the display for 20 s (either the green LED (READY) or the red LED (WARN) will light, the LED will light for 15 s and then turn off for 5 s).
 - If you press the power on/off button during the last 5 s of the countdown, the device returns to the next countdown phase.
5. Finally, “OFF” is shown on the display, the device turns off automatically together with an audible signal.
6. For further testing, continue according to points 2–4 above.
7. Ideally, perform testing with an interval of at least 2 min.

Note: After use, turn the device off by pressing the power on/off button for 2 s.


www.sunnysoft.cz © Sunnysoft s.r.o., distributor

Measurement Errors (Warning) Notice

Error message	Cause & Solution	Error message	Cause & Solution
	If the measurement result is higher than 0.25 mg/l BrAC (0.50 ‰ per mille), “WARN” is shown and “ALARM” sounds.		The battery voltage is too low to perform a test. The device turns off automatically after showing “BATT”. Replace with 2x 1.5 V batteries.
	If the result is higher than 2.00 mg/l BrAC (4.00 ‰ per mille), “Hot” is shown together with an audible signal.		Short exhalation time during the test (“Flo”); the user did not exhale smoothly and continuously for at least 3∼5 s. Press the power on/off button to perform another test.
	If you perform more than 500 tests, “CAL” is shown. Even though the unit still works, recalibration is required. Contact the distributor/manufacturer.		If the ambient temperature is outside the operating temperature range (5∼40 ℃), “OUT” is shown and the device turns off. Keep and use the device within the recommended temperature range.

© Sunnysoft s.r.o., distributor


Battery Replacement

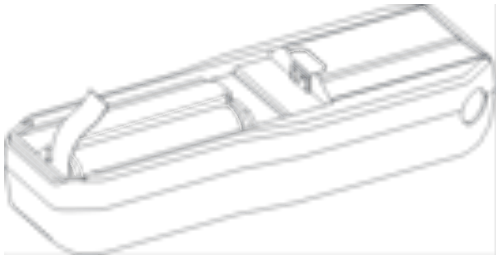
* Remove the protective tape before use.

Applying slight pressure to the battery cover located on the back of the device, push (①) and slide (②) the cover away from the device and release the cover. Insert AAA-type alkaline batteries according to the following diagram.


Step 1

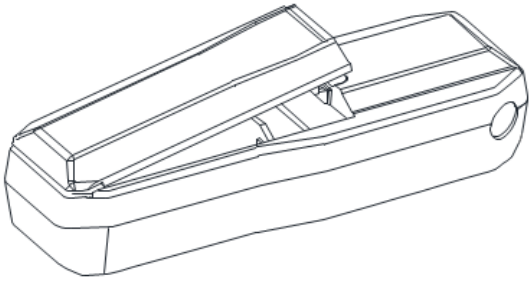
Step 2

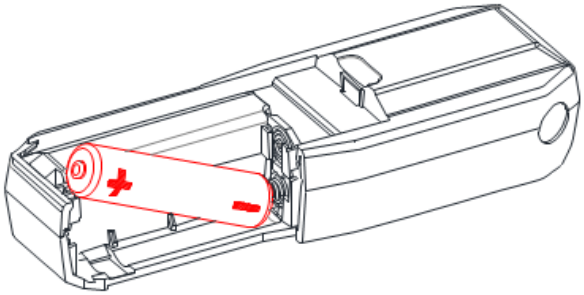
Step 3

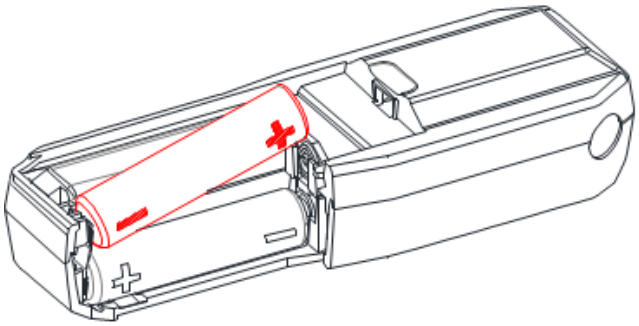Step 4

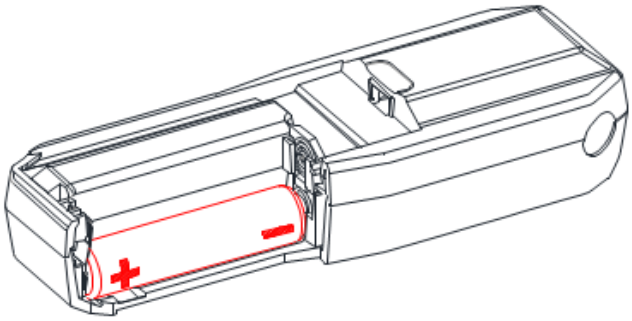Step 5

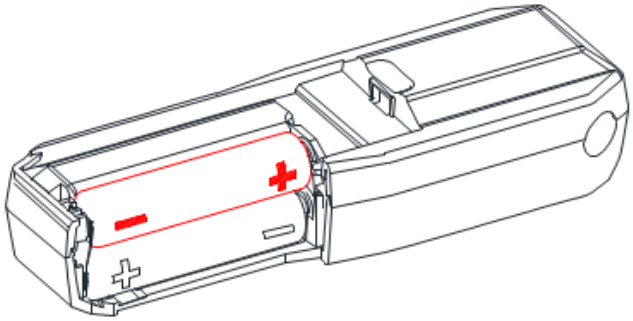Note: when inserting the batteries, observe the correct polarity. Insert from the negative terminal.

www.sunnysoft.cz © Sunnysoft s.r.o., distributor

Notice

1. Use the device approximately 20 min after eating, drinking, and smoking.
 - The reason is that absorption of alcohol from the digestive tract into the blood takes approximately 20 minutes, and residual alcohol that remains in the mouth disperses after the stated time.
2. Avoid testing in strong wind or in a closed room with polluted air.
3. If "**BAT**" appears on the display, replace the 1.5 V batteries with new ones.
4. Do not blow cigarette smoke into the device, do not pour in liquids, and do not insert food. Otherwise there is a risk of damaging the sensor.
5. Do not use the device as a tool to determine whether you or someone else should drive any motor vehicle or operate machinery. The device serves only for indicative determination.
6. Store the device at a temperature of 0–40 °C. Use the device only at a temperature in the range of 5–40 °C.
7. The correlation between the alcohol content in the breath and the alcohol concentration in the blood depends on many variables, including environmental factors (e.g. air quality, wind, humidity, temperature, etc.) and the health condition of the tested person.
8. Under normal use, the CA05FS requires recalibration every 6 months or after 500 tests to achieve accurate measurement results. Consult the distributor about recalibration.
9. Measurement results from the device cannot be used in court. The driver bears full responsibility for driving under the influence of alcohol.
10. Do not use this device in connection with an ignition interlock device or while sitting in a motor vehicle.
11. The test result may increase or decrease over time depending on many factors. Obtaining a test result at one moment does not mean that the test result will be the same a few minutes later.
12. Do not allow yourself or anyone else to drive under the influence of alcohol. If you or someone else drinks alcohol, always ask another sober driver to drive the vehicle.

www.sunnysoft.cz © Sunnysoft s.r.o., distributor

Specifications

Sensor	Advanced Fuel Cell Sensor
Accuracy	±0.005% BAC at 0.050% BAC
Detection range	0.00 ∼ 2.00 mg/l BrAC 0.00 ∼ 4.00 ‰ per mille
Warm-up time	15 ∼ 35 s
Response time	3 ∼ 22 s
Dimensions	111 mm x 35 mm x 24 mm
Weight	51.1 g (73.7 g with batteries)
Power supply	2x 1.5 V alkaline batteries
Operating temperature	5 ∼ 40 ℃
Storage temperature	0 ∼ 40 ℃
Calibration	Every 6 months or after 500 tests/uses

Distributor
Sunnysoft s.r.o.
Kovanecká 2390/1a
190 00 Praha 9
Czech Republic
www.sunnysoft.cz

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